iHerz[®]

INTELLIGENTE PFLEGE

compressor NEBULIZER SYSTEM model iFlow

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DEAR IHERZ® COMPRESSOR NEBULIZER OWNER

Thank you for choosing **iHerz®** Compressor Nebulizer model iFlow. We are sure that having appraised worthily the high quality and reliability of this device you will become a regular user of **iHerz®** products.

Your nebulizer is intended for use in the treatment of asthma, COPD and other respiratory ailments in which an aerosolized medication is required during therapy. Aerosol easily and quickly comes to both upper and lower respiratory tracts.

Before starting to use the compressor nebulizer, please, study the user manual carefully. Following the rules stated in the user manual will provide the faultless operation of the compressor nebulizer for a long period of time. For all questions concerning the unit, please, contact a service centre of **iHerz®** distributor in your country.

TABLE OF CONTENTS

INTRODUCTION	6
What Is a Nebulizer and How Does It Work?	6
What Is Nebulizing Therapy Used For?	6
What Are the Advantages of Nebulizing Therapy?	6
BEFORE STARTING TO USE THE UNIT	7
Important Safety Information	7
Working Principle	8
Unit and Accessories	9
Nebulizer Kit	10
CARRYING OUT NEBULIZING THERAPY	10
Filling the Medicine Cup	10
Operating the Unit	11
Changing the Air Filter	11
UNIT MAINTENANCE	12
Cleaning of the Mouthpiece	12
Cleaning of the Nebulizer Kit	12
Contamination by Microbes	12
Troubleshooting	13
TECHNICAL SPECIFICATIONS	14
SYMBOL DESCRIPTION	15
WARRANTY	16

INTRODUCTION

What Is a Nebulizer and How Does It Work?

The compressor nebulizer system is a device that turns a liquid drug into a fine aerosol under the influence of the compressed air. Aerosol contains large and small particles. When breathing, the gravity makes the larger particles accumulate on the walls of the upper respiratory passages while the smaller particles have enough time to reach the lower bronchial branches. Therefore the drug comes to all parts of the patient's bronchial tree, including the smallest bronchi and alveoli, and gets into the blood vessels on the walls of the alveoli in big therapeutic doses with tiny or none side effects.

What Is Nebulizing Therapy Used For?

Nebulizer systems are used by the patients with:

- acute respiratory diseases,
- obstructive pulmonary disease,
- chronic inflammatory processes in respiratory passages, such as chronic bronchitis, bronchial asthma (incl. chronic or unstable ones, exacerbation of asthma), chronic pharyngitis,
- respiratory passage infections,
- abnormalities in secretion quality and secreting in bronchi,
- different types of coughing.

What Are the Advantages of Nebulizing Therapy?

Nebulizing therapy gives a possibility to quickly act upon the inflammatory zone by big drug doses which increases the influence of the treatment and reduces the risk of side effects to other organs. Nebulizing therapy is the easiest way of inhalation treatment and it is also an economical means of treatments from the point of view of purposeful use of drug. Aerosol treatment with the help of nebulizing systems can be used for children, adults and infirm patients.

BEFORE STARTING TO USE THE UNIT

Important Safety Information

- Carefully read the instructions before use. Keep this instruction manual for future references.
- Use the device only as an inhalator for therapeutic purposes. Medicine type, mode and time of application should be prescribed by your doctor. Installation operations have to be done according to the instruction in this handbook and only for the uses declared by the constructor.
- This device is unsuitable for use with inflammable anaesthetic mixtures containing air, oxygen or nitrogen protoxide.
- Always disconnect the device from the main socket after use and before cleaning operations.
- Do not touch the device if it accidentally falls into the water. Unplug it immediately and do not use afterwards. Contact one of authorized service centers.
- Do not touch the device with wet hands.
- During use, make sure that the unit is placed on a flat and stable surface.
- Do not use the device when a patient is sleeping.
- Contact the service center for any maintenance, do not open the device. If the device does not work properly turn it off and consult the user manual.
- Do not leave any packaging materials (plastic bags, cardboard boxes, and so on) within reach of children to prevent dangers.
- Keep accessories out of reach of children. Children and people with reduced physical or mental capacities need to use the device only under close surveillance of an adult who has read this instruction manual. Keep the medicine cup out of reach of children under 36 months because it contains small parts that could be ingested.
- Always disconnect the device from the main socket when it is not used.
- It is recommended to use one personal medication set of accessories.
- Do not leave pharmaceutical residual in the medicine cup and/or in the accessories at the end of the therapy.
- Cleaning and/or maintenance must be done only after turning off the device and unplugging the power cord from the outlet.
- Do not pull the cable or appliance in order to remove the plug.
- This is a medical device for use on medical prescription. Must be run as indicated in this user instructions. It's important that patients read and understand the information for use and maintenance of the unit. Contact

your dealer or service centre for any questions.

- If the equipment's plug does not fit the socket, contact a qualified technician to substitute the plug. Do not use multiple sockets and/or extensions.
- Do not obstruct the air intake.
- The installation must be done according to the instructions and use declared by the manufacturer. Improper installation and/or improper use may cause damage to people, animals or things. The manufacturer can not be held responsible.
- Do not use the device while bathing and showering.
- Do not expose the device to the action of dust, heat (that could deform the housings), sunlight or other atmospheric agents that could compromise the safety of the device itself.
- Do not use the device nearby sources of potential interference (example: heat sources like fireplaces, electric or gas heaters; steam sources like boilers or kettles, electromagnetic radiation sources like microwave ovens, Wi-Fi devices).
- Keep away from pets that could contaminate the surface in contact with the patient or in general hinder the right performance of the therapy.
- The appliance generates compressed air. Side effects could come from the intolerance to the medicine used.

Electric devices must be disposed of in accordance with the locally applicable regulations, not with domestic waste.


KEEP THIS INSTRUCTION MANUAL DURING THE WHOLE SERVICE LIFE OF THE UNIT

WARNING: Please use only those drugs that were prescribed to you by your doctor, **DO NOT** execute self-treatment and **DO NOT** change the methods/frequency prescribed without consulting your doctor.


Working Principle

This unit is compact, simple to operate and convenient to carry. It uses an oil-free, single-cylinder piston pump to produce air pressure. The airflow created by the compressor comes out of the jet outlet and then enters the


nebulizer kit. The airflow pathway into and out of the compressor unit is illustrated in the figure below.

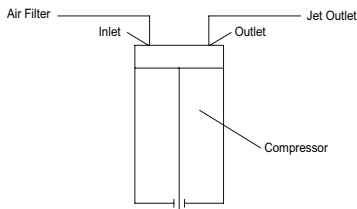Unit and Accessories

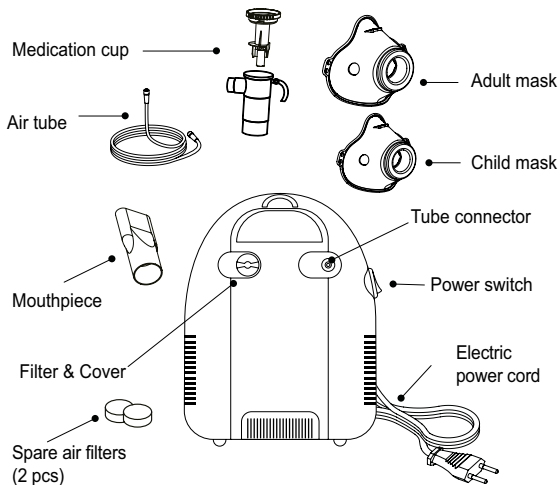Nebulizer Kit

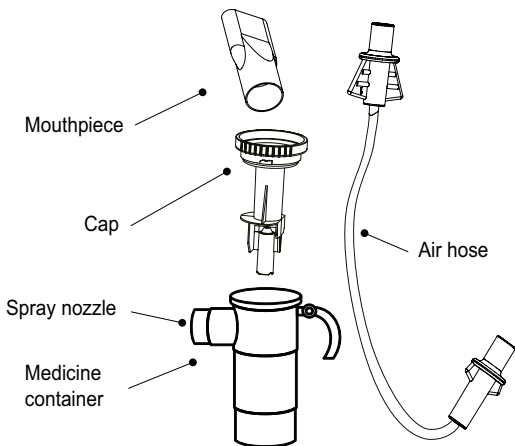
CARRYING OUT NEBULIZING THERAPY

Filling the Medicine Cup

Follow the cleaning instructions in this manual under "Unit Maintenance" prior to using your nebulizer for the first time or after it has been stored for an extended period of time.


1. Place the compressor on a flat, stable surface within reach.
2. Gently open the medication cup unscrewing its cover (figure A).
3. Add the prescribed amount of medication into the medication cup. Be sure that the medicine does not exceed the upper notch of 8 mL (figure B).
4. Reassemble the nebulizer cup by carefully twisting the medication cup

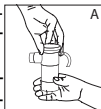and cover together. Make sure that the two parts fit securely (figure C). The stem inside the medication cup inserts into the cap's neck.

5. Connect one end of the air tube to the medication cup and the other end to the connector on the device (figure D).

6. Install the needed accessory on the medicine cup: adult or pediatric mask, or mouthpiece (figure E). Note: The mouthpiece is to be used for oral inhalation only.

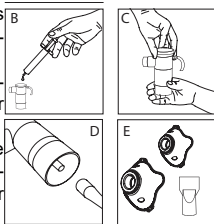
Operating the Unit

To operate the Nebulizer, do the following:

- Plug in the compressor but do not connect the nebulizer kit yet. Make sure that the air filter is attached to the compressor.
- Turn on the power switch.
- After the compressor runs properly for 1-2 minutes, connect the jet outlet of the compressor to the inlet of the nebulizer kit by the air hose.
- Begin treatment according to your physician's instructions. Inhale and exhale naturally. Be sure to stop when the medicine is used up.
- When you have finished, turn the compressor unit off before you disconnect the air hose.

Changing the Air Filter

The filter should be replaced after 100 hours of use or when it turns grey. For filter replacing, carefully pull the lid out of its housing. Remove the old filter and fit in a new one. Install the filter holder back in its housing, making sure that it is securely fixed. Use the filters supplied by the manufacture only.

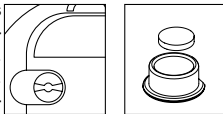RECOMMENDATIONS FOR CARRYING OUT NEBULIZING THERAPY:

- Sit down in a comfortable and upright position.
- If a mask (adult or pediatric) is used, be sure that it securely adheres to the face. Mask has holes for the air flow during the expiration.

- **Breathe easy and have a short pause at the end of each breath.**
- **Therapeutic application ends when the medicine comes out in a discontinuous mode (spluttering).**
- **When using the mouthpiece, do not tilt the medicine cup over 45°.**

UNIT MAINTENANCE

After each therapy remove the device from the socket and clean the device and the accessories. Attention: use only the accessories supplied with the product, other accessories can lead to ineffective therapy.

Cleaning of the Mouthpiece

Use drinking water to clean carefully this accessory. It can be sterilized by immersion in a chemical disinfectant (oxidizing electrolytic chloride), following disinfectant manufacturer's instructions and rinsed off with drinking water. Leave to dry naturally in a clean place.

Cleaning of the Nebulizer Kit

After use, open two parts, wash with soap and hot water or soak in warm water for 15 minutes and rinse it carefully. Remove excess of water with a soft, clean and dry cloth and leave to dry naturally in a clean place. For more thorough, hygienic cleaning, add some vinegar into the water. It is recommended to use one personal medicine cup for each patient, to be used for 6 months or 120 treatments. Do not use boiling water to clean the nebulizer kit. Carefully wash masks and hoses with disinfectant solution only.

Contamination by Microbes

In case of pathologies of microbial or viral origin, where infection is possible, the device and its accessories should not be shared among patients.

Troubleshooting

PROBLEM	FACTOR	SOLUTION
Too much noise	No filter in the unit.	Put a filter in the unit.
The nebulizer does not run	Compressor switch is off.	Turn the unit on.
	Nebulizer was not cleaned after its previous use.	Clean the unit.
	Air hose is bent.	Remove any bends or kinks in the air hose.
	The filter is clogged.	Replace the filter.
	No medication is left.	Add the appropriate amount of medication prescribed by your physician to the cup.
	The cap is not attached properly.	Reattach the cap.
Water drops have formed in the air hose.	The medication was overfilled.	Add appropriate amount of medication to cup, connect the air hose to the air compressor, and turn the unit on.
	The nebulizer kit was washed but not dried.	Open and close the nebulizer kit repeatedly until the drops disappear.

NOTE: If the device is not functioning correctly in spite of inspections effected, contact the iHerz® service center in your country.

TECHNICAL SPECIFICATIONS

Model	iFlow
Power	230 V AC, 50/60Hz
Power consumption	below 65W
Nebulizer flow	6 L
Compressor Flow	≥8 L
Medicine container capacity	2-8 mL
Particle size	0.5 to 6 microns
MMAD	2.44 microns
Average nebulization rate	≥0.2 mL/min
Residue amount	0.15 mL
Noise level	below 65dBA
Operation mode	30 minutes on, 30 minutes off
Max pressure range	≥32 Psi, ≥220 Kpa, ≥2.2 Bar
Operation pressure range	≥10 Psi, ≥70 Kpa, ≥0.7 Bar
Operating temperature range	10 °C to 40 °C (50 °F to 104 °F)
Operation humidity	10% to 95% RH
Storage temperature range	-20 °C to 70 °C (-4 °F to 158 °F)
Storage humidity range	10% to 95% RH
Dimension (L x W x H)	167 x 142x 93 mm
Weight	1.46 Kg
Protection against electric shock	Class II equipment
Degree of safety in the presence of flammable anesthetics or oxygen	No AP/APG (not suitable for use in the presence of flammable anesthetics or oxygen)
Contents	Compressor nebulizer machine, medication cup, adult mask, pediatric mask, mouthpiece, air tube, spare air filters (2 pcs)

This device has an expected service life of about 2.000 therapy cycles, based on an average therapy duration of 8 min. Electromedical devices require special care. It is essential that the devices are installed and used in accordance with the manufacturer's specifications, taking into account the requirements of electromagnetic compatibility. There is a potential risk of electromagnetic interference with other devices, in particular with other devices for analysis and treatment. Radio and mobile telecommunications devices or portable RF equipment (mobile phones or wireless connections) may interfere with the operation of electromedical device.

The specifications and the outer design are subject to changing for improvement reasons without previous notice.

SYMBOL DESCRIPTION

SYMBOL	MEANING	SYMBOL	MEANING
	Alternating current		CE mark
	Double insulation		Do not use the appliance when taking a bath or shower
	BF applied part	ON	Device is ON
	Obligation to read instruction	OFF	Device is OFF
	General warning symbol		Keep dry
IP 21	Protection against solid and liquid bodies penetration (protected against solid bodies over 12 mm; protected against access with a finger; protected against vertically falling water drops).		Temperature limitation

	This product (including its related accessories) must be delivered at the end of its service life to an electronic waste disposal/recycling facility by the user or dealer.		Manufacturing date
	Manufacturer		European representative
	Serial number		Direct current

WARRANTY

The guarantee of this device is valid for 1 year from the date of purchase. The warranty does not cover consumable materials, and those parts subject to normal wear and tear. The warranty is only valid upon presentation of the warranty card with purchase date and stamp of trade organization and receipt.

Should the unit be damaged due to user misuse or negligence, the manufacturer or authorized dealers shall not be responsible for any loss or damage to the unit. All assembly, extension, adjustment or repair of this unit must be performed only by the authorized service center of **iHerz®** distributor in your country. For repairs or spare parts purchase, please, contact the authorized service center of **iHerz®** distributor in your country.

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